

MEDLINK *SDR+*

Digital X-Ray Made Simpler



MEDLINK IMAGING


a **VIEWWORKS** company

Features with Patient Comfort in Mind

The Medlink SDR+, is a digital radiography (DR) system designed for a wide range of general and specialized diagnostic imaging.

The Semi-automated, pre-programmed positioning allows for maximum efficiency, allowing for up to 60 customized, pre-programmed positions of any exam type.

The enclosed Digital flat panel detector offers flexible imaging of all anatomy with superior resolution at industry leading lower doses. The Medlink SDR+ Integrated 10" touch screen display console automatically controls all system components while providing the operator with unprecedented simplicity and ease-of-use. Whether working in a general radiology or orthopedic application, the Medlink SDR+ system provides the user with customized workflow, unmatched functionality and exceptional diagnostic quality.

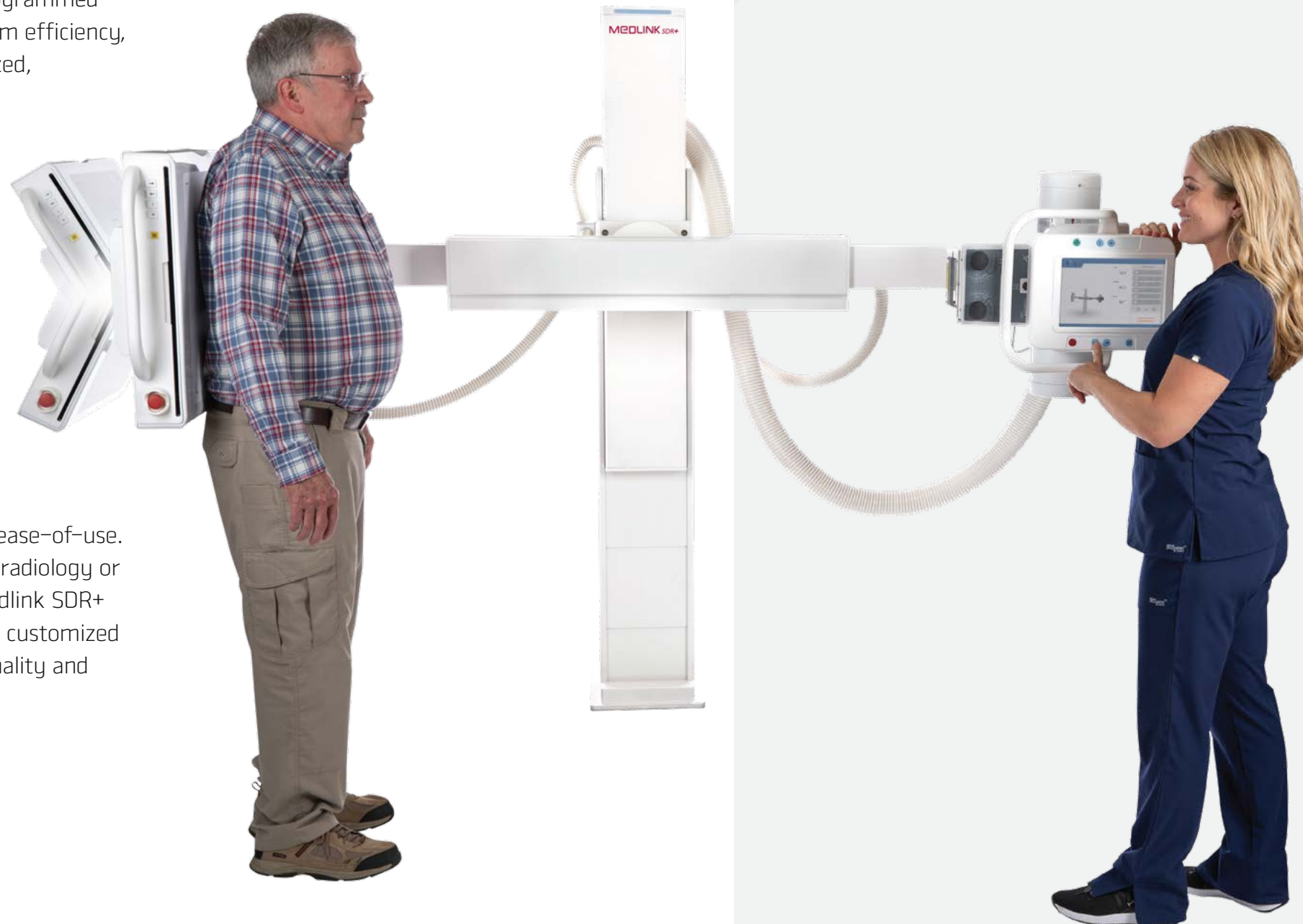

MEDLINK SDR+

Offers the following features:

- Intuitive 10" touchscreen display simplified system operation
- Up to 60 semi-automatic pre-set positions of stand height, orientation, and SID
- Removable grid
- Integrated High Frequency X-ray Generator, with customized Anatomical Programming (APR) from the digital technologist workstation
- Integrated high performance digital flat panel DR detector, with fast image display and exceptional image quality
- Occupies the smallest footprint within an 8" ceiling without limiting full-functionality
- Technologist friendly features with collimator and system movement and collision sensors
- Integrated safety features with multiple emergency-stop buttons and collision sensors
- Optional: Lightweight mobile patient imaging table, (max. Patient weight 516lbs.)
- Optional (AEC) Automatic Exposure Control
- Dual laser LED collimator

The Acuity SDR+ is also offered in our Standard configuration. Offering the same versatility of the Acuity SDR+, but with our standard non-touchscreen display and advanced features. The Acuity SDR Standard offers the same exceptional digital imaging quality as the SDR+ but with standard features at a budget friendly price point.

Simplified Operation at your Fingertips with our Intuitive 10" Touchscreen



Unlimited Exam Possibilities

Quick, easy stand positioning along with motorized source-image-distance (SID) adjustment enables a wide spectrum of upright imaging of the chest, abdomen and other skeletal exams as well as supine imaging of patients on the optional mobile table including various oblique techniques. Ergonomic user interface with multi-language support, simple APR technique selection with automated exposure parameter adjustment settings for patient size, enable a high patient throughput. The Medlink SDR+ offers a unique feature-rich design making it universally useful in all applications of diagnostic radiography.



Workstations

Featuring the VIVIX-S 1717V DR solution and a ready-to-go plug and play Dell Monitor.

Workstations includes 5-year manufacturer's warranty.

VXvue Viewing Software

The VXvue acquisition software is an award-winning intelligent software solution which utilizes the most advanced imaging processing. VXvue's image viewer program contains specialized modules to acquire, adjust, and manage diagnostic images created by the VIVIX-S Series.



The Medlink SDR System

Standard Features, Budget-Friendly Price

The Medlink SDR is also offered in our Standard configuration. Offering the same versatility and exception digital imaging quality as the Medlink SDR+, with our standard non-touchscreen display, no programmable auto-positions and limited workflow features.



Medlink SDR

MEDLINK IMAGING

a viewworks company

Linking Technology to the Imaging Professional

For the past 25 years, Medlink Imaging, LLC has been servicing the medical imaging market as a master distributor of radiographic products. Medlink has been leading the charge to expand from individual panel sales to full-room solutions. We are able to offer some of the highest-quality, competitive equipment packages available in today's market. From analog products to complete DR solutions, Medlink is the one-stop destination for all your radiographic needs.

MEDLINK I M A G I N G

a viewworks company