

# ClearRay™ **2000**

Integrated with MasterX™ System




**THE ONLY 4-WAY ELEVATING TABLE ON THE MARKET**




**7" monitor at Tube head  
for image review**



**Linak elevating column**

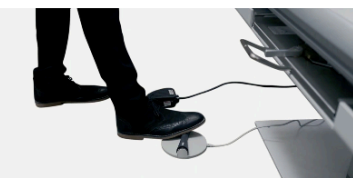

**Foot controlled**

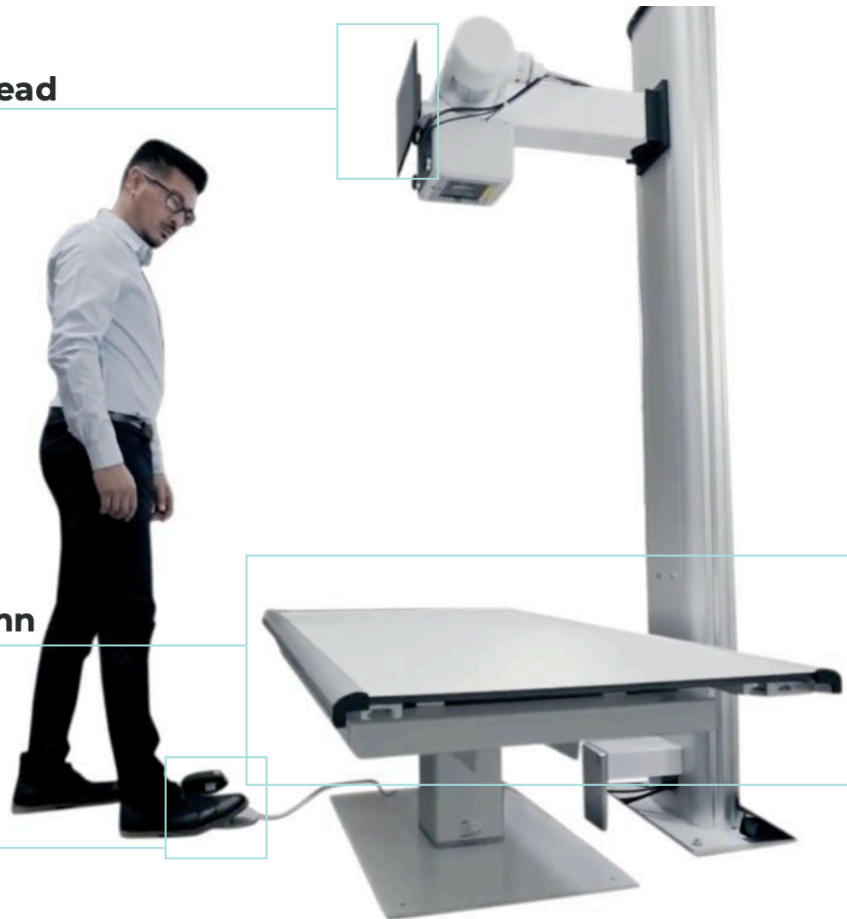

# The only 4-way elevating table on the market.

Open table platform for convenience and easy cleaning.

Sturdy high quality construction with tube column able to fit low ceilings.

One touch foot control function for easier patient positioning.

Table elevation from 18" to 31" from floor.

Electromagnetic Breaking System of table.

Single LAN cable connection — No power supply needed.

# Voyance®

By Medicatech USA®

An industry leading complete software suite that delivers a secure and cost-effective way to acquire, view and share DICOM Images. Voyance® is one of the easiest, most user friendly interfaces. Its modern design & attention to the user experience along with its features & functions make it a great tool to every user.

**Instant installation and auto-configuration**

**Easy, user-friendly Modern user interface**

**Interface consistency across acquisition and viewing platforms**

**Optimum Workflow designed for efficiency**

**ClairVoyance® Image Processing Algorithm**

**Supports a wide range of detectors and generators**

**X-RAY Generator integration (CPI & DRGEM)**

**Auto Stitching and Auto Cropping**

**Side by side comparisons**

**Fully DICOM 3.0 Compatible:  
MWL, MPPS, and all DICOM ops.**

**Reject Analysis**



# MasterX™ | Tethered

CASSETTE SIZE PREMIUM MULTI-FRAME 17"X 17" DETECTOR

POWERED BY VOYANCE®

User-friendly & ergonomic design  
Minimum dose, maximum diagnostic accuracy  
Extraordinary patient throughput & fast workflow



Single LAN Cable



Dual Energy  
Subtraction



Fast Workflow



FOCUS Scintillator  
technology

# Maximum Diagnostic Accuracy

## Fast

image preview  
within 1 second.


## Optimized

for x-ray system  
integration


## Improved

patient throughput

Others

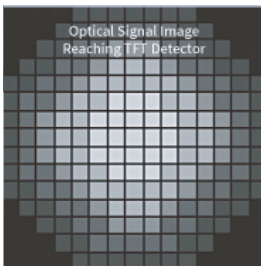

FOCUS  
Technology

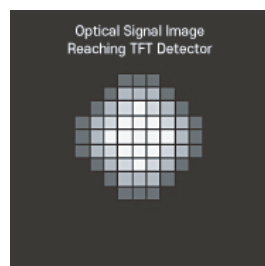


## FOCUS Scintillator Technology

FOCUS (Fluorescent Optical Csl Upgraded Structure) has more focused optical signal compared to conventional Csl. As a result, ultra low noise is achieved by FOCUS technology.

## Single LAN Cable

With a single Lan cable, installation and maintenance of detectors become more simplified. The elimination of interface box allows easier and faster installation.

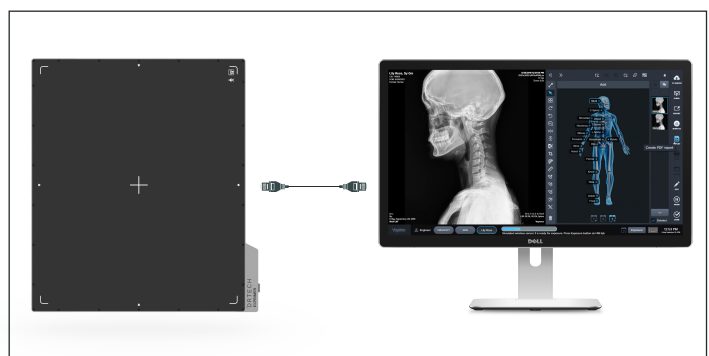

DRTECH  
EXPRIMER